




Y-Knotless™ Flex

All-Suture Anchor

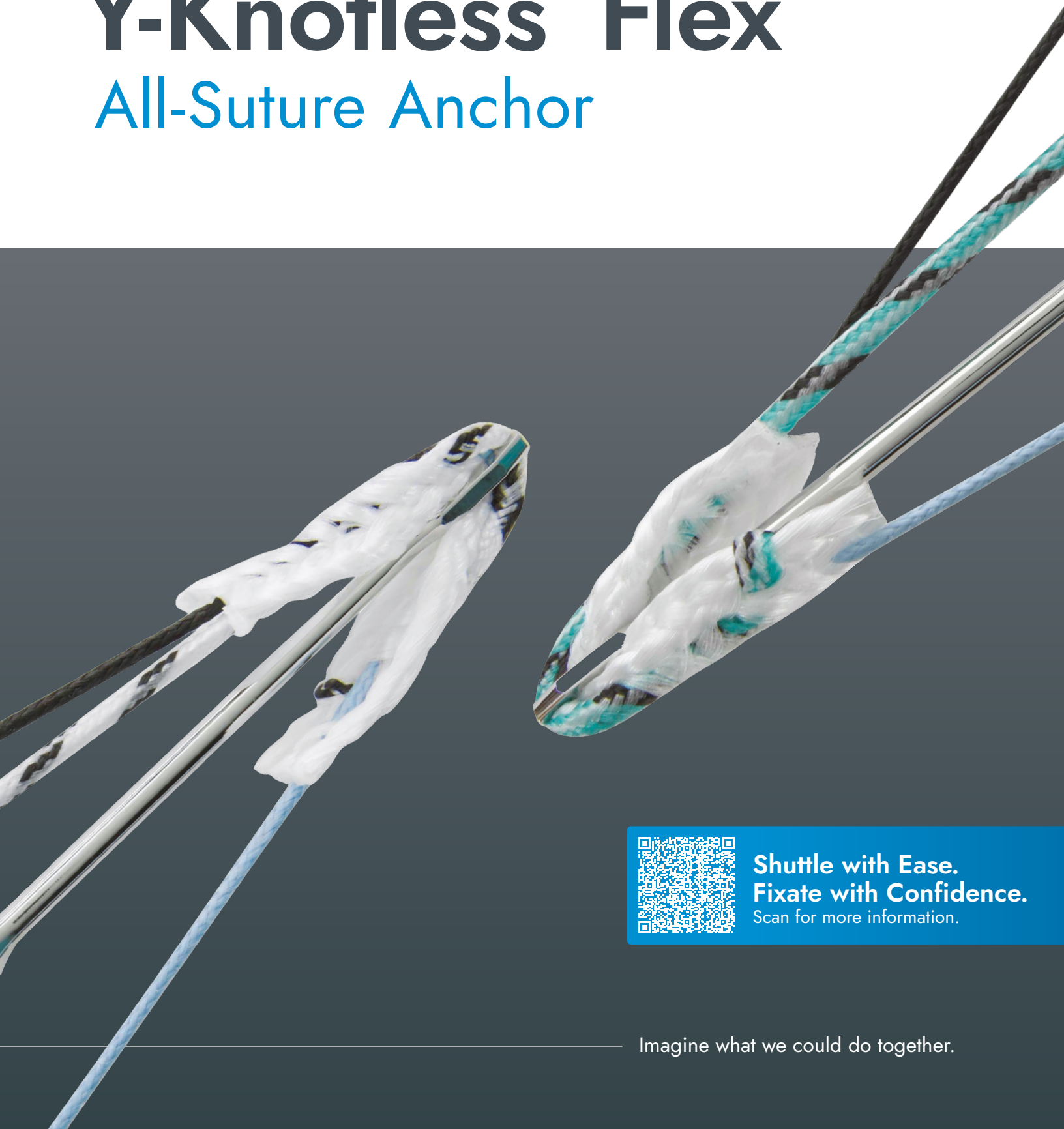

**Shuttle with Ease.
Fixate with Confidence.**
Scan for more information.

Imagine what we could do together.



Y-Knotless™ Flex All-Suture Anchor for Labral Repair in the Glenoid

With next generation knotless technology, Y-Knotless™ Flex all-suture anchor enables controlled labral tensioning with a 1.8 mm drill, helping to minimize bone removal while maximizing fixation strength. Built on the trusted 360° FormFit™ radial expansion platform, the Y-Knotless™ Flex anchor delivers consistent fixation while streamlining surgical workflow.

Designed for Confidence:

- Engineered for effortless shuttling, requiring 46% less force to shuttle the knotless mechanism compared to the leading competitor¹
- Dual-color shuttling suture provides unmistakable visual clarity throughout the procedure
- Two high-contrast repair suture configurations designed to streamline multianchor labral repair
- Proven 360° FormFit™ radial expansion enables secure fixation in a 1.8 mm pilot hole



Imagine what we could do together.


Engineered for Effortless Shuttling

The Y-Knotless™ Flex all-suture anchor features a thickened suture segment that supports controlled opening of the knotless mechanism during shuttling, and a gradual taper in the repair suture to facilitate smooth passage through the anchor. Together, these design elements reduce resistance during deployment and support efficient shuttling.

46% reduction in shuttling force¹

when mechanically tested compared to the leading competitor. This design supports a predictable and consistent shuttling experience.

Peak Shuttling Force (N)

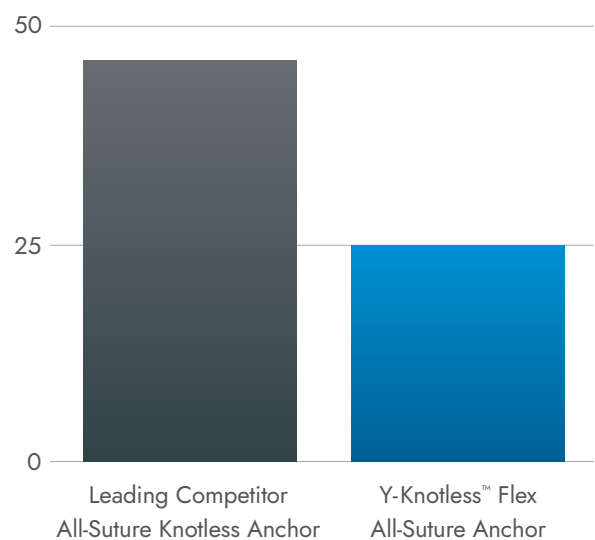

Imagine what we could do together.

Color-Coded Clarity at Every Step


Dual-color shuttling suture transitions from blue to black, enabling immediate visual identification of the looped and pull strands both arthroscopically and in the operative field.

Striped-to-solid foldover indicator enables consistent, reproducible foldover. Even in a wet environment.

Two fixation suture color options simplify every procedure.

Striped fixation suture pattern provides clear visual indicator of suture movement during tensioning.

#2 Hi-Fi® braided, non-absorbable repair strand engineered for strength and reduced tissue tear-through, demonstrating over 50% lower cutting rate versus FiberWire.²

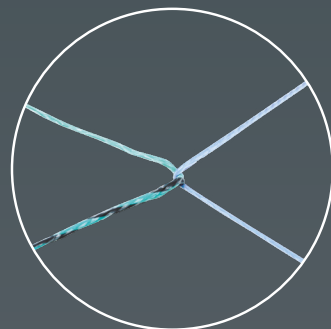


Radially Secure for Stronger Repairs

The Y-Knotless™ Flex all-suture anchor leverages proven 360° FormFit™ radial expansion to create subcortical fixation with characteristics similar to traditional solid-body anchors, while maintaining the small footprint of an all-suture design.




Fixation Strength¹ (N)

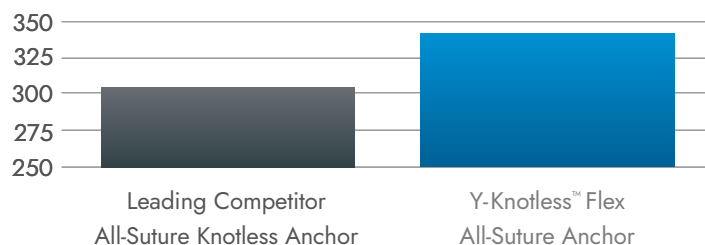


Reinforced Forked-Tip Inserters

A thoughtfully engineered forked-tip anchor inserter features reinforced construction designed to reduce the risk of in-bone breakage during anchor placement.

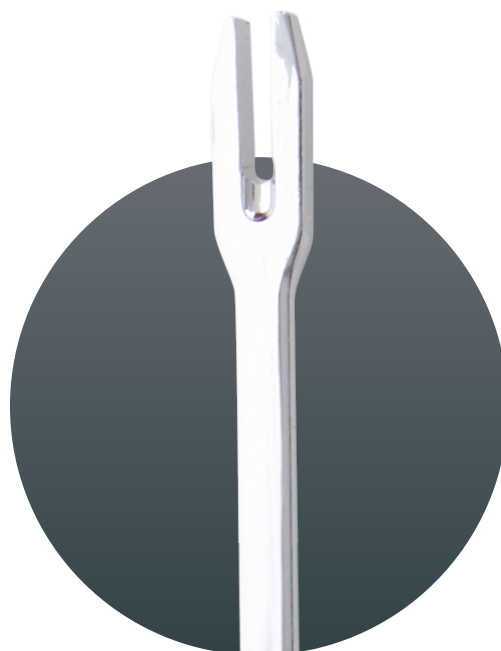

Imagine what we could do together.


Surgical Technique

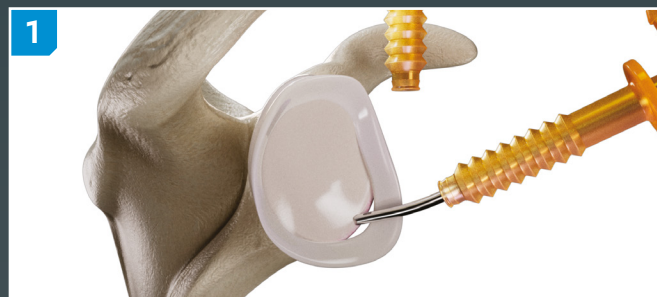


1 Position drill guide on the glenoid rim and drill a 1.8 mm pilot hole until the drill collar contacts the guide handle, providing a positive stop.




2 Insert the Y-Knotless™ Flex all-suture anchor through the drill guide and impact until the inserter handle is flush with the guide.

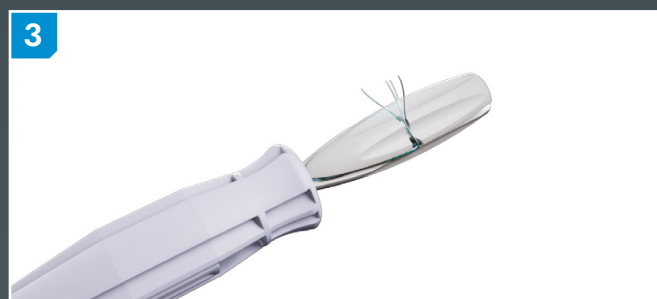


3 Uncleat sutures and remove anchor inserter.

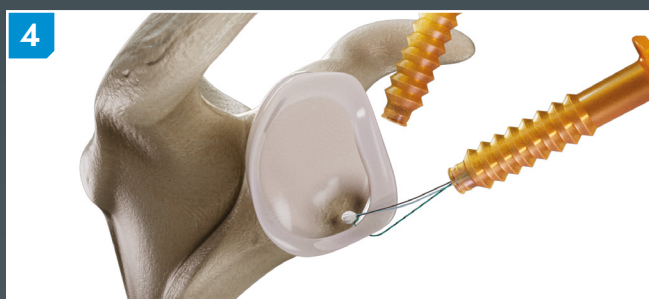


4 Set the anchor by pulling all three suture tails to achieve subcortical fixation via 360° FormFit™ radial expansion.

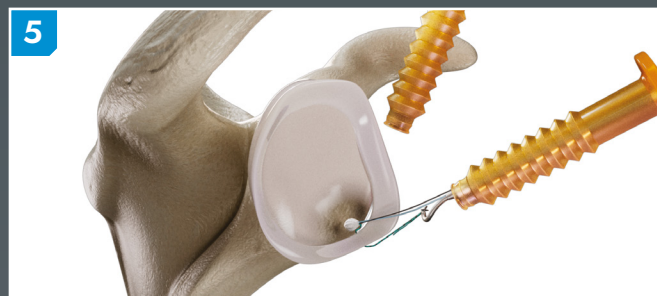

5 Retrieve the striped repair suture through the labrum.

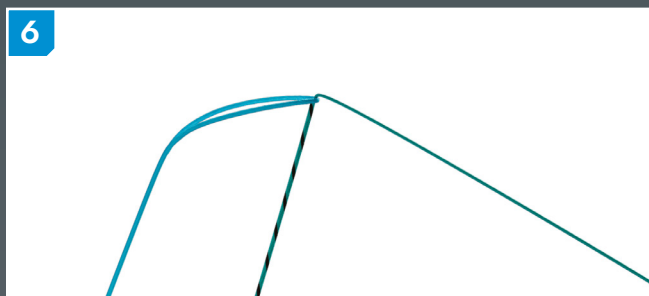

6 Pass the repair suture through the blue loop and fold at the foldover indicator (striped-to-solid transition).

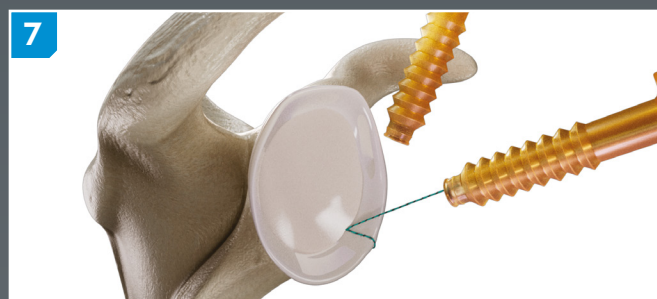

7 Shuttle the repair suture by pulling the black suture limb until the repair suture passes through the locking mechanism and exits the cannula.



8 Tension the repair suture to the desired level and repeat as needed for circumferential labral repair.

Imagine what we could do together.

Ordering Information



Implants	Catalog Number
Y-Knotless™ Flex Anchor, Green/Black #2 Hi-Fi® Suture, 1.8 mm	YKF-G
Y-Knotless™ Flex Anchor, White/Black #2 Hi-Fi® Suture, 1.8 mm	YKF-W

Imagine what we could do together.

Ordering Information



Disposable Instrumentation	Catalog Number
1.8 mm Y-Knot® Disposable Instrument Kit with Curved Guide Y-Knot® Disposable Curved Guide Y-Knot® Disposable Obturator, Blunt Y-Knot® 1.8 mm Disposable Drill Bit	YKF-DGC
1.8 mm Y-Knot® Disposable Instrument Kit with Straight Guide Y-Knot® Disposable Straight Guide Y-Knot® Disposable Obturator, Blunt Y-Knot® 1.8 mm Disposable Drill Bit	YKF-DGS



Imagine what we could do together.

Ordering Information

Reusable Instrumentation	Catalog Number
1.8 mm Disposable Drill Bit	Y18D
1.8 mm Rigid Hard Bone Drill Bit	Y18RDHB
1.8 mm Drill Guide, Fishmouth	Y-G005
1.8 mm Drill Guide, Crown	Y-G006
1.8 mm Curved Guide	Y-CG01



Drill Guide, Fishmouth
Y-G005



Drill Guide, Crown
Y-G006



Drill Guide, Curved
Y-CG01

Access Instrumentation	Catalog Number
Blunt Obturator	Y-G003
Sharp Trocar	Y-G004
Disposable Flexible Obturator, Blunt	Y-OBT1
Disposable Flexible Trocar, Sharp	Y-OBT2
1.8 mm Disposable Percutaneous Pack w/ T-Guide Distension Needle Nitinol Wire 14ga Cannulated Switching Stick Cannulated Stylet 1.8 mm Disposable Guide w/ Obturator	Y-PERC18
Disposable Percutaneous Pack Distension Needle D Nitinol Wire	Y-PERC

1. CONMED Corporation. Data on file (PDD2929454). Largo, FL; 2025

2. Deranlot J, Maurel N, Diop A, Nourissat G, et al. Abrasive Properties of Braided Polyblend Sutures in Cuff Tendon Repair: An in Vitro Biomechanical Study Exploring Regular and Tape Sutures. Arthroscopy 2014; 30:1569-73.



Imagine what we could do together.

CONMED Corporation | 11311 Concept Blvd., Largo, FL 33773 | www.CONMED.com



Service & Support